



TECHNISCHE  
UNIVERSITÄT  
WIEN  
Vienna University of Technology

VIENNA UNIVERSITY OF TECHNOLOGY

FACULTY OF PHYSICS

LABORATORY III

---

# Laboratory Report

Electron Spin Resonance

---

**Authors:**

Raul Wagner

Martin Kronberger

**Group 301**

**Supervisor:**

Someone

conducted on:  
04 June 2025

# 1 Resonanceabsorbption of a passive HF-Oscillator

## 1.1 Setup

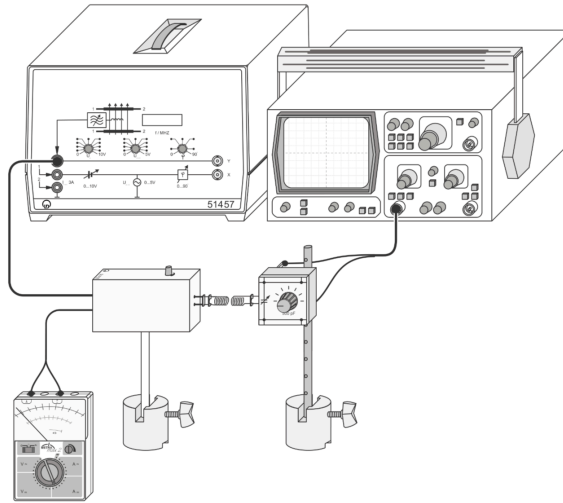


Figure 1: Experimental setup

- Connect the ESR base unit to the ESR operating unit via a 6-pin cable and set the rotary potentiometer to maximum sensitivity.
- Plug in the 30–75 MHz plug-in coil and connect the ammeter to output I via an adapter cable (measurement range 100  $\mu\text{A}$ ).
- Position the coil of the passive resonant circuit coaxially opposite the plug-in coil and connect via a BNC/4 mm measurement cable to channel I of the dual-channel oscilloscope.

## 1.2 Procedure

- Set the variable capacitor of the passive resonant circuit to position  $\text{Skt.} = 3/6$ .
- Adjust the minimum frequency on the ESR base unit.
- At the operating frequency, measure and record:
  - the frequency,
  - the voltage  $U_2$  of the “passive” coil on the oscilloscope,
  - and the voltage  $U_1 = 56\text{k}\Omega \cdot I_1$  of the RF coil.
- Increase the frequency stepwise and repeat the measurement.
- Perform additional measurement series with  $\text{Skt.} = 2/6$  and  $1/6$ .
- Remove the passive resonant circuit and record another measurement series.

**1.3 Measurement values**

freq / MHz	$U_2$ / V	$U_1$ / V
22.1	1.17	2.48
22.3	1.25	2.31
22.5	1.3	1.9
22.7	1.29	1.75
22.9	1.25	1.59
23.1	1.11	2.34
23.3	0.93	2.67
23.5	0.76	2.86
24	0.48	3.1
24.5	0.32	3.25
25	0.25	3.34
25.5	0.19	3.37
26	0.15	3.43
27	0.09	3.54
28	0.07	3.67
29	0.05	3.75
30	0.03	3.83

Table 1: Tab. 1: Spannungen  $U_2$  und  $U_1$  bei Skt. = 3/6

# Laboratory Work III - Electron Spin Resonance

freq / MHz	$U_2$ / V	$U_1$ / V
22.5	0.49	2.97
23	0.55	3.01
23.5	0.64	3.04
24	0.79	3.02
24.2	0.87	2.97
24.4	0.96	2.84
24.6	1.08	2.72
24.8	1.15	2.36
25	1.2	1.89
25.2	1.19	1.49
25.5	1.15	1.53
25.7	1.1	2.01
25.9	1.02	2.37
26.1	0.88	2.72
26.3	0.76	2.93
26.5	0.64	3.08
27	0.45	3.33
27.5	0.34	3.43
28	0.27	3.54
29	0.17	3.67
30	0.12	3.72
31	0.09	3.82
32	0.06	3.91
33	0.05	3.96
34	0.04	4.04
35	0.04	4.05

Table 2: Tab. 2: Spannungen  $U_2$  und  $U_1$  bei Skt. = 2/6

Laboratory Work III - Electron Spin Resonance


freq / MHz	$U_2$ / V	$U_1$ / V
23	0.28	3.1
24	0.28	3.23
25	0.28	3.32
26	0.28	3.42
27	0.29	3.52
28	0.31	3.6
29	0.34	3.68
30	0.39	3.71
31	0.48	3.7
31.5	0.54	3.7
32	0.64	3.6
32.5	0.77	3.32
32.7	0.83	3.1
32.9	0.88	2.79
33.1	0.91	2.36
33.3	0.91	1.79
33.5	0.89	1.44
34	0.86	1.47
34.2	0.85	2.04
34.4	0.8	2.59
34.6	0.76	2.82
34.8	0.67	3.14
35	0.6	3.36
35.5	0.46	3.7
36	0.36	3.86
36.5	0.3	3.99
37.5	0.22	4.14
38.5	0.15	4.26
39.5	0.12	4.3

Table 3: Tab. 3: Spannungen  $U_2$  und  $U_1$  bei Skt. = 1/6

freq / MHz	$U_1$ / V
22.2	3.14
25	3.47
30	3.92
35	4.26
40	4.48

Table 4: Tab. 4: Spannung  $U_1$  ohne Resonanzabsorption

### Laboratory Work III - Electron Spin Resonance

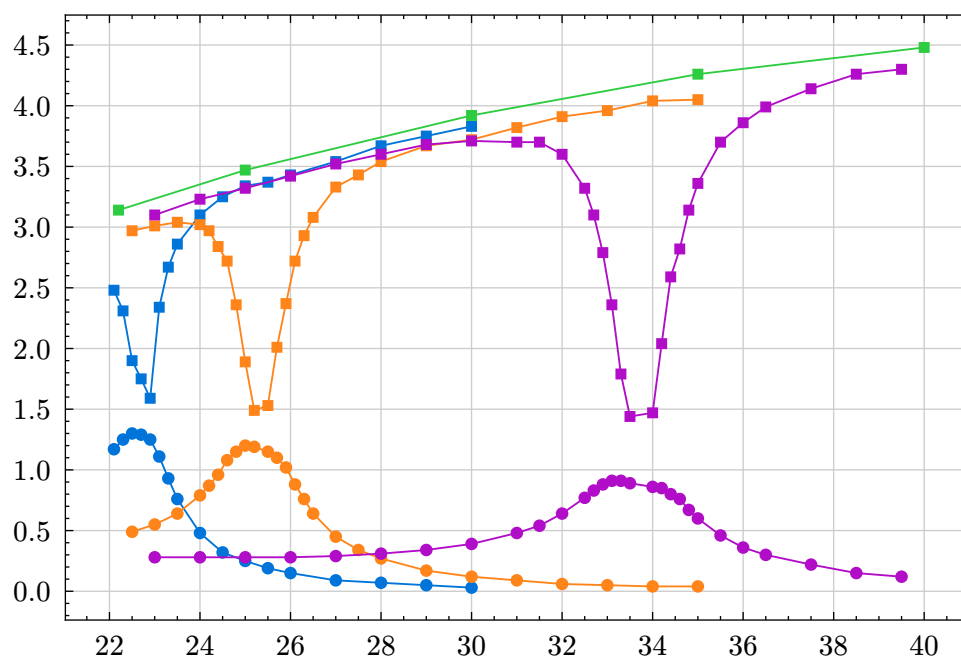

#### 1.4 Data

## 2 Electrons spin resonance on DPPH

### 2.1 Setup

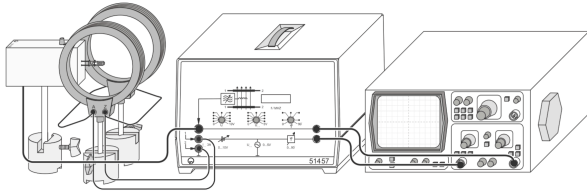


Figure 3: Experimental setup

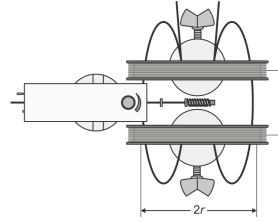


Figure 4: Experimental setup

The experimental setup is shown in Fig. 4 and 5.

- Place the Helmholtz coils parallel to each other at a center distance of 6.8 cm (equal to the mean radius  $r$ ).
- Connect both Helmholtz coils in series with the ammeter to the ESR operating unit.
- Connect the ESR base unit to the ESR operating unit via a 6-pin cable.
- Connect output Y of the ESR operating unit via a BNC cable to channel I of the dual-channel oscilloscope, and output X to channel II.

### 2.2 Procedure

#### Determination of the Resonance Magnetic Field $B_0$

- Insert the 15–30 MHz plug-in coil and place the DPPH sample centrally.
- Switch on the ESR base unit and position it so that the plug-in coil with DPPH sample is in the center of the Helmholtz-coil pair (see Fig. 5).
- Set the resonance frequency  $\nu = 15$  MHz.
- Set the modulation amplitude  $U_{\text{mod}}$  to the second scale division.
- Set the phase shift to  $0^\circ$ .
- Operate the oscilloscope in dual-channel mode:
  - Dual on
  - Time base  $2 \frac{\text{ms}}{\text{cm}}$
  - Amplitude I and II  $0.5 \frac{\text{V}}{\text{cm}}$  AC
- Slowly increase the DC voltage  $U_0$  to the Helmholtz coils until the resonance signals are equidistant (see Fig. 6).
- Switch the oscilloscope to XY mode and adjust the phase shift so that the two resonance peaks coincide (see Fig. 3).
- Vary  $U_0$  until the resonance signal is symmetric, keeping the modulation voltage as low as possible.
- Measure the DC current  $2I_0$  through the Helmholtz-coil pair and record it together with the resonance frequency  $\nu$ .
- Increase  $\nu$  by 5 MHz and adjust  $U_0$  to reestablish resonance.
- Again measure and record the current  $2I_0$ .
- Continue raising  $\nu$  in 5 MHz steps (switch to the 30–75 MHz coil at 30 MHz, and to the 75–130 MHz coil at 75 MHz) and repeat the measurements.

#### Determination of the Half-Width $\delta B_0$

- Operate the oscilloscope in XY mode:
  - Amplitude II  $0.5 \frac{\text{V}}{\text{cm}}$  AC
- Reestablish the resonance condition for  $\nu = 50$  MHz (middle plug-in coil).

- Vary the modulation voltage  $U_{\text{mod}}$  until the resonance trace spans the full screen width (10 cm) in the X-direction.
- Switch the ammeter to AC mode and measure the effective current  $2I_{\text{mod}}$  corresponding to  $U_{\text{mod}}$ .
- Increase the X-deflection, read off the width  $\Delta U$  of the resonance peak at half its height, and record it.

### 2.3 Measurement values

$\nu$ / MHz	$2 I_0$ / A	Steckspule
15	0.26	klein
20	0.35	klein
25	0.44	klein
30	0.51	klein
30	0.51	mittel
35	0.6	mittel
40	0.69	mittel
45	0.77	mittel
50	0.86	mittel
55	0.94	mittel
60	1	mittel
65	1.1	mittel
70	1.2	mittel
75	1.3	mittel
75	1.3	groß
80	1.35	groß
85	1.45	groß
90	1.55	groß
95	1.6	groß
100	1.7	groß
105	1.8	groß
110	1.9	groß
115	1.95	groß
120	2	groß
125	2.1	groß
130	2.2	groß

Table 5: Tab. 1: Stromstärke  $2I_0$  in Abhängigkeit von der Frequenz  $\nu$  des Wechselfeldes



$\nu$ / MHz	$B_0$ / mT
15	0.55
20	0.74
25	0.93
30	1.08
30	1.08
35	1.27
40	1.46
45	1.63
50	1.82
55	1.99
60	2.12
65	2.33
70	2.54
75	2.75
75	2.75
80	2.86
85	3.07
90	3.28
95	3.38
100	3.6
105	3.81
110	4.02
115	4.12
120	4.23
125	4.44
130	4.65

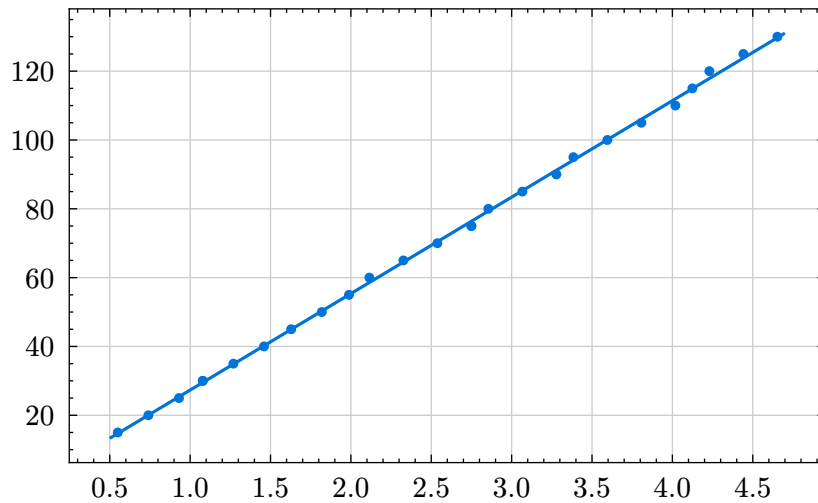
Table 6: Tab. 2: Magnetfeld  $B_0$  in Abhängigkeit von der Frequenz  $\nu$  des Wechselfeldes

Figure 5: Some Caption

## Laboratory Work III - Electron Spin Resonance

Slope:  $28.03 \frac{\text{MHz}}{\text{mT}}$

Gfactor:

$$g = 2.0026$$

from literature:  $g = 2,0036$

### **2.4 Data**