



TECHNISCHE
UNIVERSITÄT
WIEN
Vienna University of Technology

VIENNA UNIVERSITY OF TECHNOLOGY

FACULTY OF PHYSICS

LABORATORY III

Laboratory Report

Electron Spin Resonance

Authors:

Raul Wagner

Martin Kronberger

Group 301

Supervisor:

Someone

conducted on:
04 June 2025

1 Resonanceabsorbption of a passive HF-Oscillator

1.1 Setup

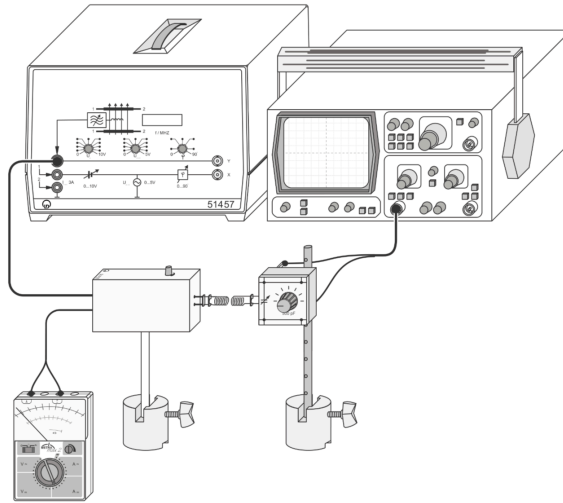


Figure 1: Experimental setup

- Connect the ESR base unit to the ESR operating unit via a 6-pin cable and set the rotary potentiometer to maximum sensitivity.
- Plug in the 30–75 MHz plug-in coil and connect the ammeter to output I via an adapter cable (measurement range 100 μA).
- Position the coil of the passive resonant circuit coaxially opposite the plug-in coil and connect via a BNC/4 mm measurement cable to channel I of the dual-channel oscilloscope.

1.2 Procedure

- Set the variable capacitor of the passive resonant circuit to position $\text{Skt.} = 3/6$.
- Adjust the minimum frequency on the ESR base unit.
- At the operating frequency, measure and record:
 - the frequency,
 - the voltage U_2 of the “passive” coil on the oscilloscope,
 - and the voltage $U_1 = 56\text{k}\Omega \cdot I_1$ of the RF coil.
- Increase the frequency stepwise and repeat the measurement.
- Perform additional measurement series with $\text{Skt.} = 2/6$ and $1/6$.
- Remove the passive resonant circuit and record another measurement series.

1.3 Measurement values

| freq / MHz | U_2 / V |
|------------|-----------|
| 11.5 | 1 |
| 12 | 1.01 |
| 12.5 | 1.15 |
| 13 | 1.2 |
| 13.5 | 1.4 |
| 14 | 1.6 |
| 14.5 | 1.8 |
| 15 | 2.2 |
| 15.5 | 2.35 |
| 16 | 2.2 |
| 16.5 | 2 |
| 17 | 1.8 |
| 17.5 | 1.25 |
| 18 | 1 |
| 18.5 | 0.8 |
| 19 | 0.7 |

Table 1: Tab. 1: Spannungen U_2 und U_1 bei Skt. = 3/6


| freq / MHz | U_2 / V |
|------------|-----------|
| 11.5 | 0.8 |
| 12.5 | 0.85 |
| 13.5 | 0.9 |
| 14.5 | 0.97 |
| 15.5 | 1.05 |
| 16.5 | 1.2 |
| 17 | 1.3 |
| 17.5 | 1.4 |
| 18 | 1.5 |
| 18.5 | 1.6 |
| 19 | 1.7 |
| 19.5 | 1.6 |
| 20 | 1.5 |
| 20.5 | 1.35 |
| 21 | 0.95 |
| 21.5 | 0.8 |

Table 2: Tab. 2: Spannungen U_2 und U_1 bei Skt. = 2/6

| freq / MHz | U_2 / V |
|------------|-----------|
| 18 | 0.75 |
| 19 | 0.8 |
| 20 | 0.85 |
| 21 | 0.9 |
| 22 | 0.95 |
| 23 | 1 |
| 23.5 | 1.05 |
| 24 | 1.1 |
| 24.5 | 1.15 |
| 25 | 1.1 |
| 25.5 | 1.05 |
| 26 | 1 |
| 26.5 | 0.95 |
| 27 | 0.9 |
| 27.5 | 0.8 |

Table 3: Tab. 3: Spannungen U_2 und U_1 bei Skt. = 1/6

Es konnte ebenso die Messreihe mit keinem passiven nicht durchgeführt werden!

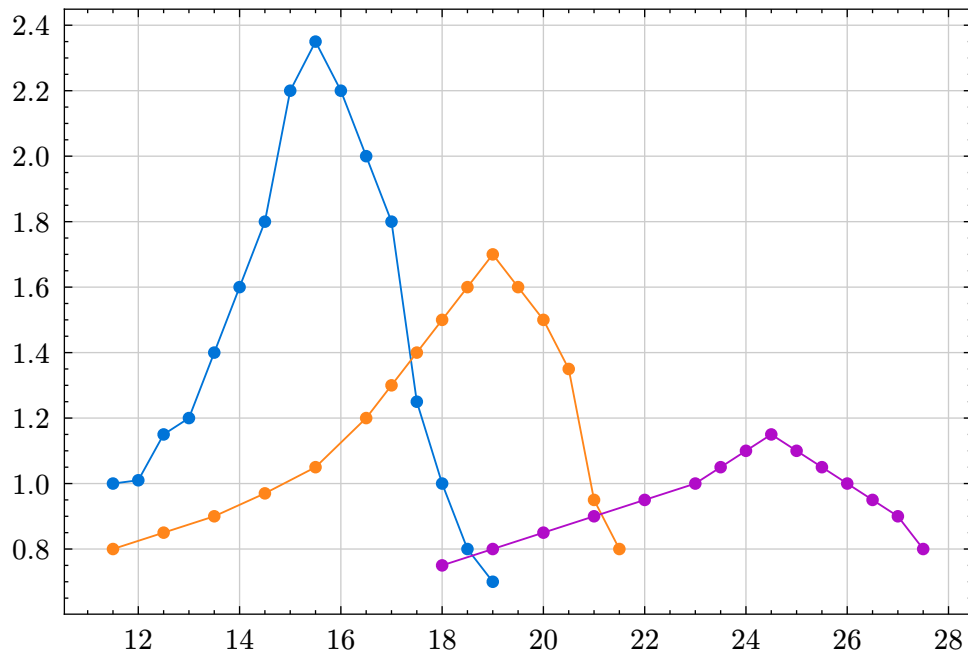

1.4 Data

2 Electrons spin resonance on DPPH

2.1 Setup

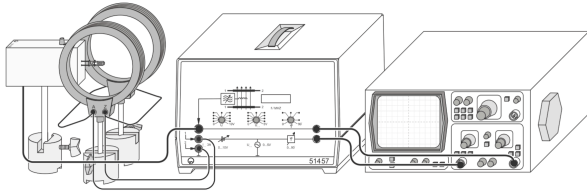


Figure 3: Experimental setup

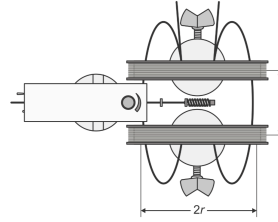


Figure 4: Experimental setup

The experimental setup is shown in Fig. 4 and 5.

- Place the Helmholtz coils parallel to each other at a center distance of 6.8 cm (equal to the mean radius r).
- Connect both Helmholtz coils in series with the ammeter to the ESR operating unit.
- Connect the ESR base unit to the ESR operating unit via a 6-pin cable.
- Connect output Y of the ESR operating unit via a BNC cable to channel I of the dual-channel oscilloscope, and output X to channel II.

2.2 Procedure

Determination of the Resonance Magnetic Field B_0

- Insert the 15–30 MHz plug-in coil and place the DPPH sample centrally.
- Switch on the ESR base unit and position it so that the plug-in coil with DPPH sample is in the center of the Helmholtz-coil pair (see Fig. 5).
- Set the resonance frequency $\nu = 15$ MHz.
- Set the modulation amplitude U_{mod} to the second scale division.
- Set the phase shift to 0° .
- Operate the oscilloscope in dual-channel mode:
 - Dual on
 - Time base $2 \frac{\text{ms}}{\text{cm}}$
 - Amplitude I and II $0.5 \frac{\text{V}}{\text{cm}}$ AC
- Slowly increase the DC voltage U_0 to the Helmholtz coils until the resonance signals are equidistant (see Fig. 6).
- Switch the oscilloscope to XY mode and adjust the phase shift so that the two resonance peaks coincide (see Fig. 3).
- Vary U_0 until the resonance signal is symmetric, keeping the modulation voltage as low as possible.
- Measure the DC current $2I_0$ through the Helmholtz-coil pair and record it together with the resonance frequency ν .
- Increase ν by 5 MHz and adjust U_0 to reestablish resonance.
- Again measure and record the current $2I_0$.
- Continue raising ν in 5 MHz steps (switch to the 30–75 MHz coil at 30 MHz, and to the 75–130 MHz coil at 75 MHz) and repeat the measurements.

Determination of the Half-Width δB_0

- Operate the oscilloscope in XY mode:
 - Amplitude II $0.5 \frac{\text{V}}{\text{cm}}$ AC
- Reestablish the resonance condition for $\nu = 50$ MHz (middle plug-in coil).

- Vary the modulation voltage U_{mod} until the resonance trace spans the full screen width (10 cm) in the X-direction.
- Switch the ammeter to AC mode and measure the effective current $2I_{\text{mod}}$ corresponding to U_{mod} .
- Increase the X-deflection, read off the width ΔU of the resonance peak at half its height, and record it.

2.3 Measurement values

| ν / MHz | $2 I_0$ / A | Steckspule |
|-------------|-------------|------------|
| 30 | 0.53 | middle |
| 35 | 0.63 | middle |
| 40 | 0.71 | middle |
| 45 | 0.79 | middle |
| 50 | 0.89 | middle |
| 55 | 0.97 | middle |
| 60 | 1.06 | middle |
| 65 | 1.15 | middle |
| 70 | 1.23 | middle |
| 75 | 1.33 | middle |
| 80 | 1.41 | middle |
| 80 | 1.53 | small |
| 90 | 1.65 | small |
| 95 | 1.67 | small |
| 100 | 1.7 | small |
| 105 | 1.74 | small |
| 110 | 1.79 | small |
| 115 | 2.05 | small |
| 120 | 2.16 | small |

Table 4: Tab. 1: Stromstärke $2I_0$ in Abhängigkeit von der Frequenz ν des Wechselfeldes

Laboratory Work III - Electron Spin Resonance

| ν / MHz | B_0 / mT |
|-------------|------------|
| 30 | 1.12 |
| 35 | 1.33 |
| 40 | 1.5 |
| 45 | 1.67 |
| 50 | 1.88 |
| 55 | 2.05 |
| 60 | 2.24 |
| 65 | 2.43 |
| 70 | 2.6 |
| 75 | 2.81 |
| 80 | 2.98 |
| 80 | 3.24 |
| 90 | 3.49 |
| 95 | 3.53 |
| 100 | 3.6 |
| 105 | 3.68 |
| 110 | 3.79 |
| 115 | 4.34 |
| 120 | 4.57 |

Table 5: Tab. 2: Magnetfeld B_0 in Abhängigkeit von der Frequenz ν des Wechselfeldes

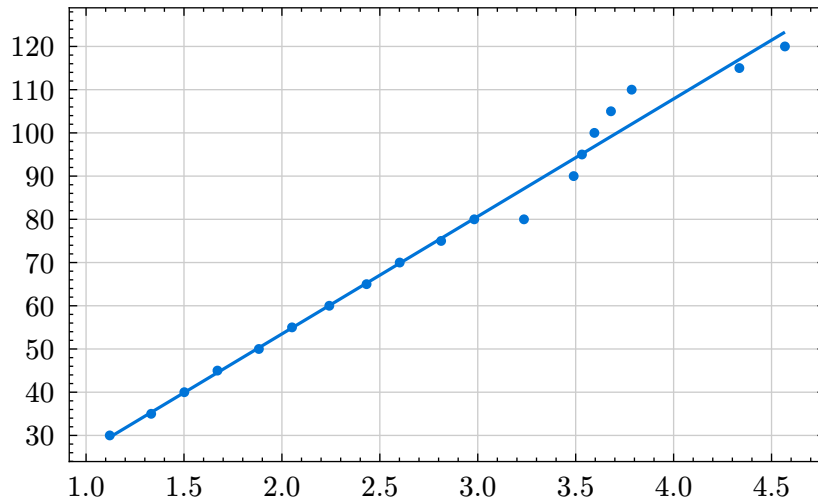


Figure 5: Some Caption

Abgelesene halbwertsbreite:

$$\delta U = 0.95V$$

$$\delta I = 0.078A$$

$$\delta B_0 = 0.33 \text{ mT}$$

$$\text{Slope: } 27.19 \frac{\text{MHz}}{\text{mT}}$$

Gfactor:

Laboratory Work III - Electron Spin Resonance

$$g = 1.9426$$

from literature:

$$g = 2.0036$$

2.4 Data